


Sterility Assurance Focus Product Guide

Healthmark offers a variety of solutions to help your department document and verify equipment and process performance. Ensure effective steam sterilization with easy-to-use and reliable sterility assurance indicators. Interested in learning about these and other solutions? Contact us!

hmark.com
800.521.6224

Sterilization Tape

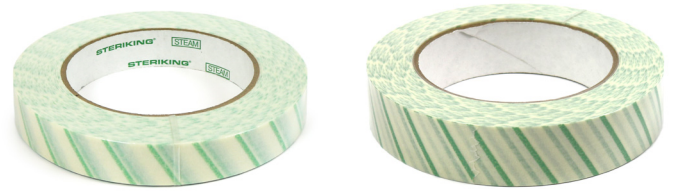

Designed for closing packaged items that require wrapping for sterilization. When correctly exposed, the Type 1 (Process Indicator) changes color from green to brown/black.

ITS19GL	Sterilization Tape - 0.75"	48 Rolls
ITS25GL	Sterilization Tape - 1.00"	36 Rolls

Mini Bowie-Dick Test Pack



This Type 2 (Bowie-Dick test indicator) is designed to detect the presence of residual air in pre-vacuum steam sterilizers operating at 132°C and 134°C for 3.5 minutes.

61301606667	Mini Bowie-Dick Test Pack	30/Box
-------------	---------------------------	--------

Bowie-Dick Test Card



This Type 2 (Bowie-Dick test indicator) is designed for daily monitoring of pre-vacuum steam sterilizers operating at 132°C for 3.5 minutes to detect the presence of residual air.

6005500396	Bowie Dick Test Card Holder	1/Box
6005500614	Bowie Dick Test Card	30/Bag

Steam Migrating Integrator



During steam sterilization, the chemical pellet of this Type 5 (Integrating Indicator) melts and migrates as a blue color along the paper wick. The migration is visible through a pass/fail window.

61301606825	Steam Migrating Integrator Strips	250/Bag
-------------	-----------------------------------	---------

ChemiSix Emulator Test Pack



When exposed to steam sterilization conditions, this Type 6 (Emulating Indicator) gives an integrated response by changing color from purple to green.

61301606235	ChemiSix Emulator Test Pack	40/Box
-------------	-----------------------------	--------